
LAPAROSCOPIC STAPLER

Trizone ™


Disposable Endoscopic Articulating Reload With FUNNEL Technology

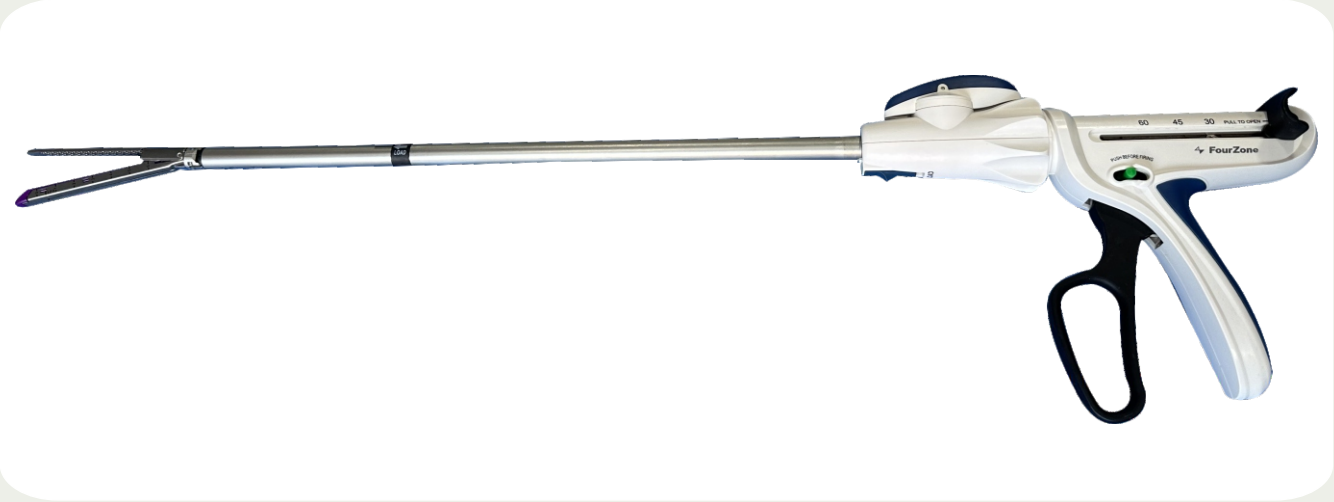

TRIZONE RELOAD AND JAW DETAIL

fourzone.in
info@fourzone.in | sales@fourzone.in

Trizone™

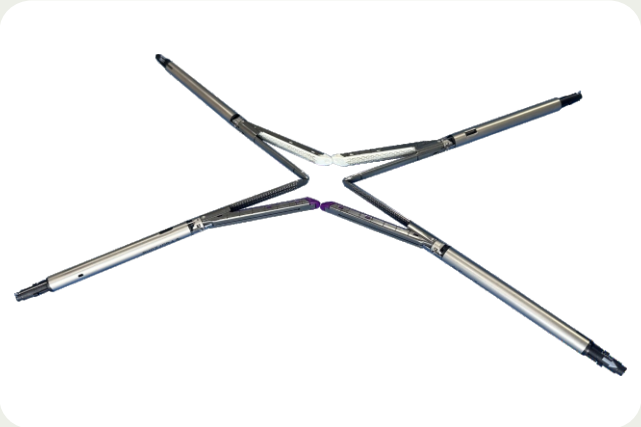



TENZR POWERED HANDLE, MANUAL STAPLER AND RELOAD CARTRIDGES

Six rows of staples are laid in three different heights, tallest at the outside. The tissue is compressed gradually across the jaw rather than uniformly, which is what keeps blood supply at the staple line without giving up hold on the thick edge.

The blade is replaced with every reload, so the first firing and the fifth are cutting with the same edge. An isometric height retainer holds the anvil geometry constant while the staples form.

Key Differentiators

01 **Six-row cross pattern**
Three staggered heights across six rows, tallest outermost.

02 **Fresh blade every reload**
The cutting edge is never reused, which keeps the transection clean and reduces tissue drag.

03 **Safety interlock**
A spent reload cannot be fired again.

04 **360° rotation, 45° swing**
Full shaft rotation with articulation both ways to reach the target site.

Ordering — stapler bodies

CODE	CARTRIDGE	STAPLE HEIGHTS (MM)	STAPLE LINE	STAPLES
FZHPL-45W	White	2.4/ 2.6 / 3.0	46 mm	66
FZHPL-45P	Purple	3.4/ 3.6 / 4.0	46 mm	66
FZHPL-45B	Black	4.4/ 4.6 / 4.8	46 mm	66
FZHPL-60W	White	2.4/ 2.6 / 3.0	61 mm	90
FZHPL-60P	Purple	3.4 /3.6 / 4.0	61 mm	90
FZHPL-60B	Black	4.4 / 4.6 / 4.8	61 mm	90

FOURZONEHEALTHCARE

Secure and stable healthcare products

FourZone Healthcare develops and manufactures devices that improve patient outcomes and advance minimally invasive surgery. Everything in this catalogue is made to one standard, whether it is used in the first minute of a case or the last.

HEADOFFICE

A-205 & 208, Premium House
Near Gandhigram Railway Station
Ashram Road, Ahmedabad 380006
Gujarat, India

info@fourzone.in
sales@fourzone.in
fourzone.in

SCAN TO VISIT